

## Determination of Solidus and Liquidus Alloys Temperatures

### Using a Differential Scanning Calorimeter

The solidus temperature is defined as the temperature at which the first deviation from the base line appears (Figure 1). The deviation signals that a phase change is taking place. Thus the solidus temperature has been reached. The tangent line, drawn in by the differential scanning calorimeter (DSC), is an approximation of the solidus temperature.

The liquidus temperature is defined as the temperature at which the graph returns to the baseline. The end of

the deviation signifies the end of the phase change, i.e., the alloy has reached the liquid phase and the liquidus temperature. The tolerance on the values from the chart is  $\pm 3^{\circ}\text{C}$ .

The appearance of the DSC output is severely affected by the scan rate. Scan rate is the speed at which the temperature is increased during the scan. A slower scan yields a "tighter" curve (i.e., it is more accurate because the sample has time to reach equilibrium).

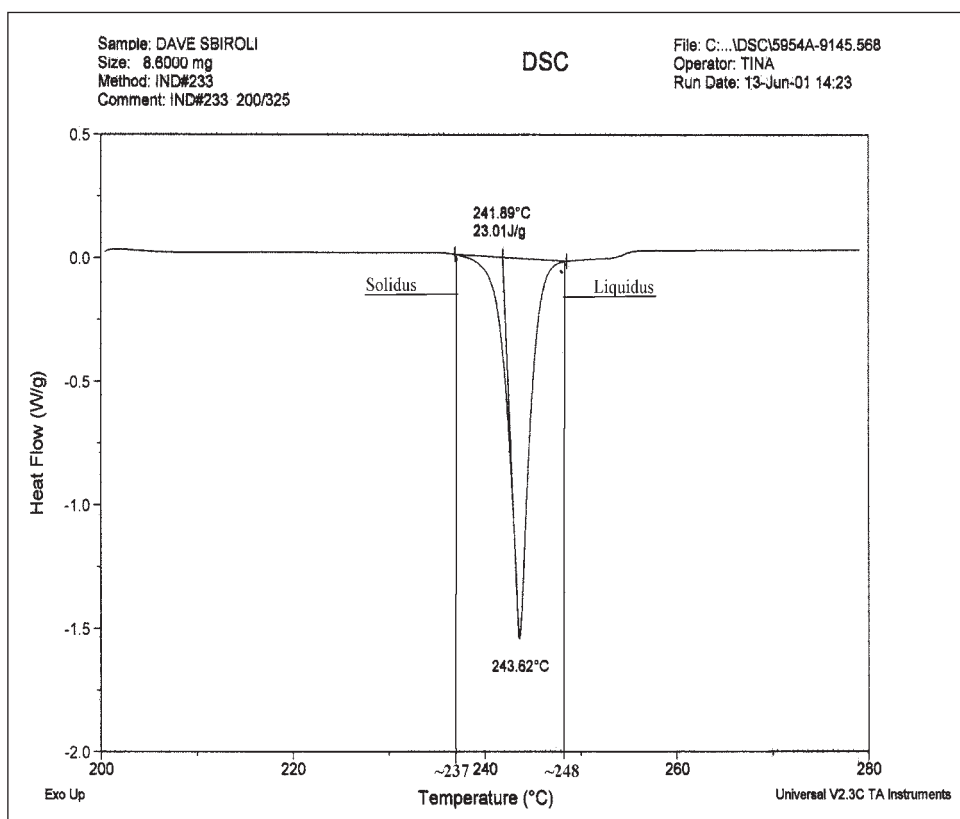


Figure 1.

This application note is provided for general information only. It is not intended, and shall not be construed, to warrant or guarantee the performance of the products described which are sold subject exclusively to written warranties and limitations thereon included in product packaging and invoices. All Indium Corporation's products and solutions are designed to be commercially available unless specifically stated otherwise.

Form No. 97831 R2

[www.indium.com](http://www.indium.com)

[askus@indium.com](mailto:askus@indium.com)

ASIA: Singapore, Cheongju, Malaysia: +65 6268 8678  
 CHINA: Suzhou, Shenzhen: +86 (0)512 628 34900  
 EUROPE: Milton Keynes, Torino: +44 (0) 1908 580400  
 USA: Utica, Clinton, Chicago, Rome: +1 315 853 4900



©2016 Indium Corporation

APPLICATION NOTE

ISO 9001  
REGISTERED